

**ACCESSORIES**




Lithotomy Crutches



Shoulder Support



Intravenous Arm board with wrist straps

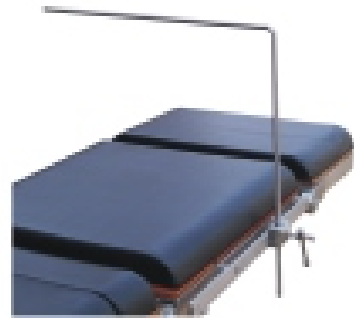

Anaesthetic Screen Holder



Foot Rest



Padded Lateral Support



Stand for Table Accessories




Swivable Mobile IV Stand



Gel/ Pu Foam Calf Support

**DRT-EH-05**

**O.T. TABLE  
3 WAY OPERATED  
ELECTRIC | BATTERY | MANUAL**

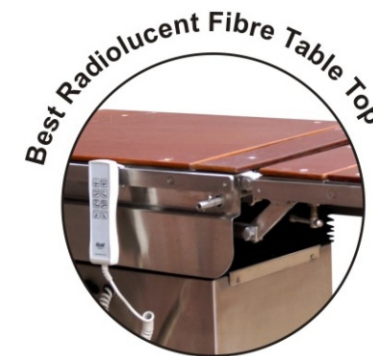

MANUAL SELECTOR



HAND HELD REMOTE AND  
OVERRIDE CONTROL PANEL FOR ALL  
FUNCTIONS



**DORNIER MEDICAL SYSTEMS**

3699, Sector 23  
Gurgaon, Haryana - 122017  
Email: [dorniermedicalsyste.ms@gmail.com](mailto:dorniermedicalsyste.ms@gmail.com)  
Ph: +91-9650047772

## TECHNICAL SPECIFICATIONS

1. Modular C-Arm image intensifier fully X-Ray radiolucent high pressure laminated scratch less electric OT table divided into 4-5 sections.
2. Table consist of advanced noiseless 100% oil free electric system for the facilities to adjust movements from hand held remote and electronic backup override control panel on column for height adjustment, trendelenburg and reverse trendelenburg, left and right tilt, back up/down and zero position for getting table back to normal position from any position.
3. The hand held green backlit remote shows symbolic presentation of all positions for clear display of operations even in dark light.
4. In the event of complete power failure, table is equipped with manual override control panel on the column for adjustment of all positions like height adjustment, trendelenburg and reverse trendelenburg, left and right tilt and back section up/down.
5. Detachable and interchangeable head and leg sections. The seat section is designed in such a way that avoids patient from slipping down even while raising the back section. Head section has facility of up & down and leg section down with split apart laterally which makes it suitable for all sizes of patients.
6. Full scratch proof table top consist detachable radiolucent, antistatic electric conductive, seamless PU molded mattress fixed through velcro at table top designed for patient comfort prevent pressure ulcer. High quality easy to clean latex free mattress of 50mm thickness.
7. Medical grade full stainless steel SS 304# thick grade base, side rails for side handles and clamps with epoxy powder coated base makes it resistant to disinfection and prevent ingress of fluid. Table material is completely shock and impact resistant.
8. Table is equipped with heavy duty antistatic four lockable corrosion proof ABS castors for easy easy swiveling and rotation of OT Table. Good base height with beveled edges allow sufficient sufficient room for surgeon leg during surgery.
9. Stand by UPS power supply rechargeable built-in maintenance free battery system with backup of 4 hours (20 operations) and can be used while charging in event of power failure.
10. The maximum weight bearing capacity of the table is 160kg.

## TECHNICAL PARAMETERS

PARAMETERS	FEATURES MODEL: DRT-EH-05
Overall Length and Width (LxW) Length without foot section	190cm (L) x 52cm (W) 137 cm
Height Adjustment(electrically)	66 cm (minimum) - 100 cm (maximum)
Lateral Tilt (electrically)	20° Left and Right
Trendelenburg & R. Trendelenburg (electrically)	25° both
Back section (electrically)	90 °up and 45° down
Leg section adjustment (manually)	30° up, 90° down and split 180°
Head section adjustment (manually)	60° up and 90° down
Optional: Built-in kidney elevation (manually)	120mm/ 5 inches
Power Source	220V - 240V, 50Hz - 60Hz

## ELECTRIC OPERATIONS



HEIGHT ADJUSTMENT  
66 TO 100cm



LATERAL TILT BOTH SIDES  
20°



TRENDELENBURG & REVERSE  
TRENDELENBURG  
25°



BACK SECTION  
+90°/ -45°